


## Connecting Options for Analytical HPLC Columns, Short Columns and Guard Columns with 2, 3, 4 und 4,6 mm ID

### 1. Guard Column and Analytical Column

- 1.1 Connection via capillary or one-piece PEEK column connector (P/N 590253)
- 1.2 Direct connection via coupling nut (P/N 424204) and stainless steel connector (P/N 424404)

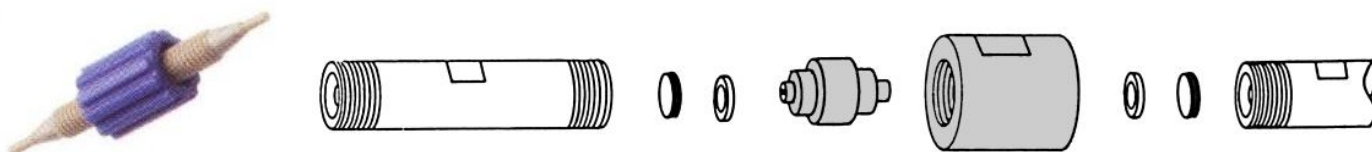


### 2. Column Connector for ECO Columns 2 - 4 mm ID-Columns

The ECO column is delivered without column heads and fitting adapters, only with end caps.

For column connection you need:

2 pcs. column head, P/N 424201

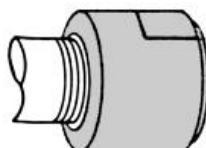


2 pcs. fitting adapter, P/N 424260



### 3. Guard Column Holder and Analytical Column

3.1 Guard column holder with UNF 10/32 female (connection via capillary tubing), P/N 424271 and 424272

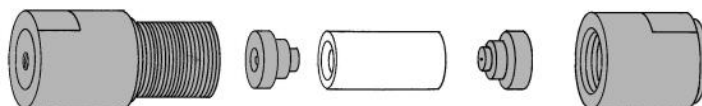

3.2 Guard column head, direct connection of guard column with analytical column, P/N 424262, 424266, and 424267

